

Medical Grade Endoscope Video Coupler




The **Medical Grade Video Coupler** is an easy to use optical component that allows one to conveniently attach any CCD Camera with the C-mount to a medical endoscope with the standard 32 mm diameter eyepiece. The objective lens of the Video Coupler relays a magnified image of the object in front of the endoscope to a CCD chip of the camera.

Video Coupler Features



FIT TO ENDOSCOPE	STANDARD 32 mm ENDOSCOPE EYEPIECE
CLAMPING	SPRING-LOADED, THREE PETAL CLAMP OR THREE BALL CLAMP
FOCUSING	ONE HAND PRECISION FOCUSING
OPTICAL WINDOWS	SCRATCH RESISTING COATING
CCD CAMERA THREAD	C-MOUNT MALE, 1.00-32
FOCAL LENGTH	25,28,30,35,40,45 mm
OPTICS	HIGH QUALITY, COATED
STERILIZATION	SOAKABLE

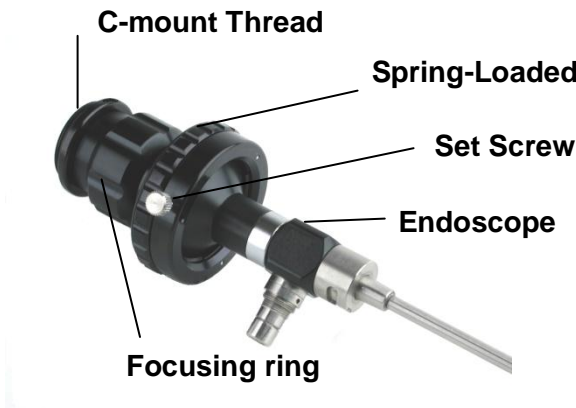

MODEL No.	Focal Length (mm)
MVC-2500	25.0
MVC-2800	28.0
MVC-3000	30.0
MVC-3500	35.0
MVC-4000	40.0



**Alternative Design
with 3-Petal Clamp**



We also design and manufacture custom video cameras, optical adapters, imaging components and systems for digital and video cameras used in medical applications.



AGS-MEDICAL, Ph: 505-5506501 or 505-5655102, Fax: 505-814-5778,

Email: sales@agstech.net, Web://www.agsmedical.com